
WATER BATH FOR SWELLING AND BOIL TESTS

SW200 & BT200

SW 200 S

BEST IN CLASS FOR:

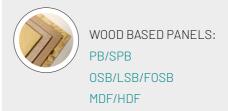

The SW200 and BT200 appliances consist of a water bath which is used to carry out laboratory and immersion tests on samples of wood-based panels (particleboard, MDF, OSB) in accordance with EN European standard requirements (it may also be used in conjunction with other IMAL equipment) such as: • Swelling and absorption (EN 317) • Boil test (EN 1087-1)

The system has a control unit for regulating temperature as required by the EN standards. Up to 2 separate timers may be set on the touchscreen display (one per sample rack). Temperature is monitored from the trend view.

The equipment comes with racks for immerging the samples for the test that is to be carried out, the 2 racks for the EN317 swelling test have been designed to house the samples in full compliance with the standard requirements and the distances indicated. The number of samples which may be immerged will depend on sample thickness plus a minimum space between one sample and the next. The SW200 can hold up to 136 samples (depending on thickness). The 2 racks for the BT200 version for the EN1087-1 test can house the samples that have been glued to the relative supports, the bath can hold 24 samples. Water temperature is set on the display in relation to the test which is to be carried out, and is regulated by a PID control; an immersion heater is


fitted inside the tank to regulate the water temperature at the required value. If the difference between water temperature set-point and the environmental temperature is less than 2°C, an optional immersion cooling circuit can be provided. Imal recommends a solution of this kind to keep the temperature where the machine is located at the recommended temperature and to ensure that the water remains "still" during the test as per the standard requirements. The system is also equipped with a level control, safety temperature control and drain tap.


The SW200 and BT200 are an integral part of the IMAL Smartlab platform. A webserver is also included to monitor the state of the unit with remote connection from PC and mobile or tablet as well.

SW BLOCK 200 and BT BLOCK 200 versions are available for carrying out tests in accordance with UIC CODE 435-2 requirements



IN CONJUNCTION WITH IB800 BOARD PROPERTY TESTER

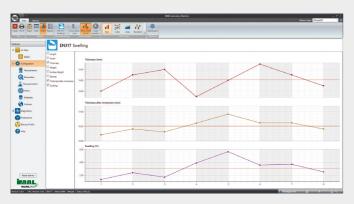

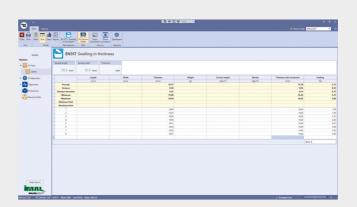

TECHNICAL DATA SW200 & BT200		
OPERATING TEMPERATURE (EN 317)	20 °C	
OPERATING TEMPERATURE (EN 1087-1)	100 °C	
TEMPERATURE STABILITY	±1°C	
RESOLUTION	0.1°C	
NOMINAL SAMPLE SIZE	50 x 50 mm	
WATER VOLUME	361	
INSTALLED POWER	2.8 kW / 230 V	
SAMPLES PER RACK SW200	24 ÷ 136 (in relation to thickness)	
SAMPLES PER RACK BT200	12	

TECHNICAL DATA SW200 BLOCK & BT200 BLOCK		
OPERATING TEMPERATURE SW200 BLOCK / BT200 BLOCK	20 °C / 100 °C	
TEMPERATURE STABILITY	±1°C	
RESOLUTION	0.1°C	
WATER VOLUME SW200 BLOCK / BT200 BLOCK	891/551	
INSTALLED POWER SW200 BLOCK / BT200 BLOCK	2.8 kW / 4.0 kW	
NUMBER OF SAMPLES PER RACK SW200 BLOCK / BT200 BLOCK	10/6	